Personvåg...3

EN Body weight scale.....7

TEKNISKA DATA

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Batterier ingår ej. Batteri: 3 x 1,5V AAA.

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SE

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Kundkontakt: +46(0)20-83 33 33 eller ICA.se

TECHNICAL DATA

Model No: EB3469. Maximum weight capacity: 180 kgs (395 lbs).

Batteries not included. Battery: 3 x 1,5V AAA.

ADVICE FOR USE & CARE

EN


1. The platform will be slippery when wet. Keep it dry!
2. Stand still while weighing.
3. DO NOT strike, shock, or drop the scale.
4. Handle the scale with care as it is a precision instrument.
5. Clean the scale with a damp cloth and avoid water entering the scale. DO NOT use chemical or abrasive cleaning agents.
6. Store the scale in a cool, dry place.
7. Always keep the scale in a horizontal position.
8. If the scale fails to turn on, check if the battery is installed or if the battery power is low. Insert or replace the battery.
9. If there is an error on the display or the scale cannot be turned off for a long time, remove the battery for about 3 seconds, then reinstall it to reset. If the problem persists, contact your dealer for assistance.
10. Do not use the scale for commercial purposes.

TECHNICAL DATA

Model No: EB3469. Maximum weight capacity: 180 kgs (395 lbs).


Batteries not included. Battery: 3 x 1,5V AAA.

EN INSTRUCTIONS FOR USE

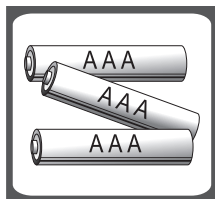BATTERY

3 X 1,5V AAA batteries (not included)

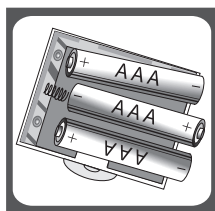TO REPLACE BATTERY

1. Open the battery case cover on the back of the scale.
2. Remove the used batteries (use a pointed object if necessary).
3. Install new batteries by inserting one side first and then pressing down the other side. Pay attention to the polarity.

OPERATION

1. To ensure accuracy, place the scale on a hard, flat surface and avoid carpets or soft surfaces.
Note: The scale is a slim version, make sure to place it on a hard surface and ensure there are no foreign objects beneath the scale before weighing.
2. Gently step on the scale to turn it on automatically. Stand evenly on the scale without moving and wait until your weight is displayed stably and locked.


TECHNICAL DATA

Model No: EB3469. Maximum weight capacity: 180 kgs (395 lbs).

Batteries not included. Battery: 3 x 1,5V AAA.


INSTRUCTIONS FOR USE

EN

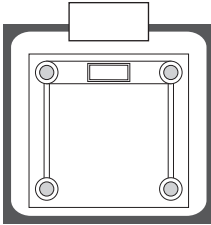AUTOMATIC SWITCH OFF

The scale will automatically switch off when you step off. Auto shut-off occurs if the display shows a stable weight reading for approximately 8 seconds.

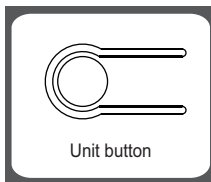TO CONVERT WEIGHT UNIT


If the scale supports both Metric (kilograms) and Imperial (pounds) systems, you can convert the weight unit as follows:

1. Locate the unit switch on the back of the scale.
2. Press the unit conversion button when the scale is powered on. The LCD will show the current weight unit. Press the button again to convert the weight unit.


LO-SYMBOL: LOW BATTERY INDICATION

Battery power is low. Please install a new battery.

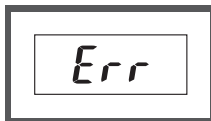ERR-SYMBOL: OVERLOAD INDICATION

The scale is overloaded (maximum weight capacity is indicated on the scale). Please step off to avoid damage.

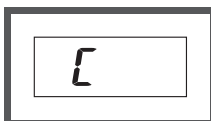C-SYMBOL: RE-TEST

If the display shows "C" when you step off the scale, it means the weight reading is incorrect. Please re-test.

TECHNICAL DATA

Model No: EB3469. Maximum weight capacity: 180 kgs (395 lbs).

Batteries not included. Battery: 3 x 1,5V AAA.

Due to the constant development of our products in terms of function and design, we reserve the right to make changes to the product without prior warning. We reserve the right to make amendments. In accordance with the WEEE directive, every member state must ensure correct collection, recovery, handling and recycling of electrical and electronic waste.

Private households in the EU can take used equipment to special recycling stations free of charge. In some member states you can, in certain cases, return the used equipment to the retailer from whom you purchased it, if you are purchasing new equipment. Contact your retailer, distributor or the municipal authorities for further information on what you should do with electrical and electronic waste.

Produced in China for ICA AB, Box 4075, 169 04 Solna, Sweden.

Customer contact: +46 (0)20-83 33 33 or ICA.se